enable lifecare
Fast intelligent care

# Configura® Comfort

## Information & Assessment Form








# USER GROUP

To provide comfortable seating to users with mild to moderate postural requirements, to help reduce the formation of pressure ulcers and help to stand or sit users with limited mobility.

## APPROPRIATE FOR:

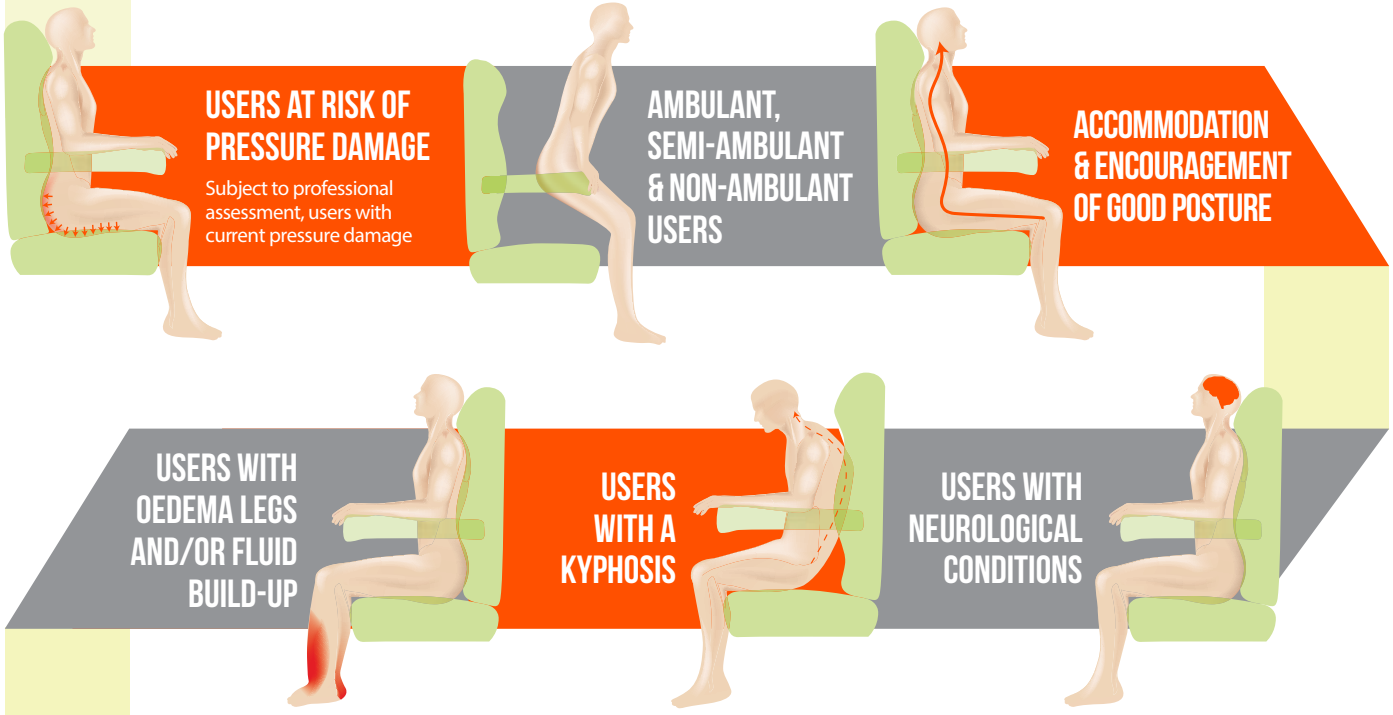


## NOT APPROPRIATE FOR:

1. Demanding postural requirements
2. Users with fixtures and contractures
3. Users with a limited extension at the knees
4. Users who need to be transported from room to room
5. Bariatric users

# INTENDED ENVIRONMENT

- Domestic use
- Care homes or nursing homes

# STANDARD CHAIR SPECIFICATION

**Configura® Comfort is the latest innovation from Accora and provides an unprecedented level of functionality and cost effectiveness.**


Comfort is the result of acting on feedback collected from our customers over the past 5 years and it includes many features you have been asking for as standard.







| |
|---|
| 'Reduce Shear' recline to maximise pelvic stability |
| Integrated self-supply overlay for easy integration of existing cushions* |
| Compatible with transfer aids including Cricket, ReTurn 7500 and Molift Raiser |
| Knock down design for easier storage and single person delivery |
| Pressure reducing fabric on backrest and armrest as standard |
| Colour matched vapour permeable fabric |
| Castors for easy manoeuvrability |
| More depth adjustment options using backrerst brackets. Reduces the need for depth adjustment kits |
| 60kg legrest capacity |

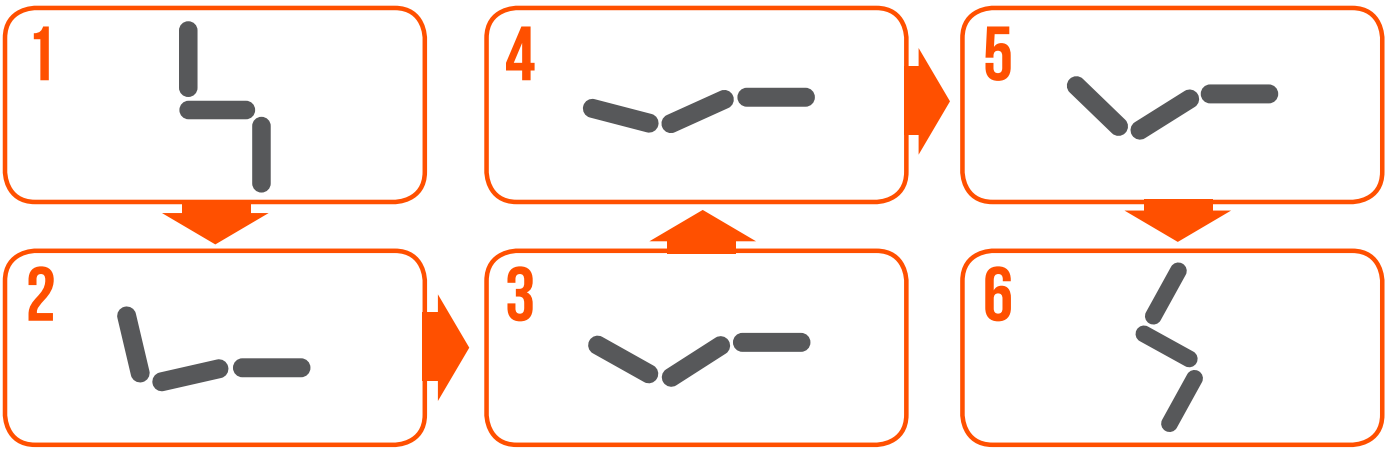


**For further information on the integrated self-supply overlay please see page 11.*

# PRESCRIBER DETAILS

Name: .....

Department: .....

Address: .....

.....

.....

Telephone: .....

Email: .....

# USER DETAILS

Name: .....

Address: .....

.....

.....

Telephone: .....

Email: .....

User weight: .....

Pressure Care: .....

Posture: .....

Mobility: .....

Diagnosis: .....

Prognosis: .....

Transfer method: .....

Current seating: .....


Problems with current seating: .....

Notes: .....

.....

.....


.....

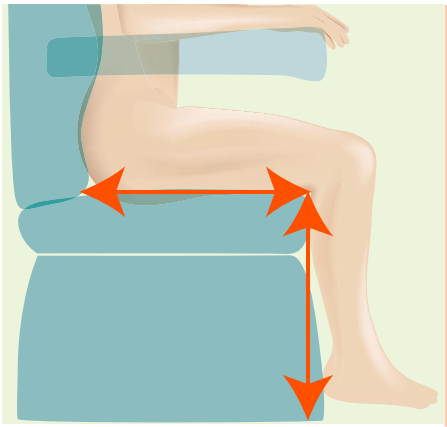

# SIZE


**Key thoughts:**

- 3 measurements are required for selecting a chair; seat depth, seat height and seat width – see below drawings for clarity.
- When selecting a chair size, prioritise the chair measurements in the following order: seat depth, seat height then seat width. This will ensure the selected size is the right overall proportion for the user
- If you cannot find the size you need, check out the other chairs in our range or contact us as we may be able to supply a special size.

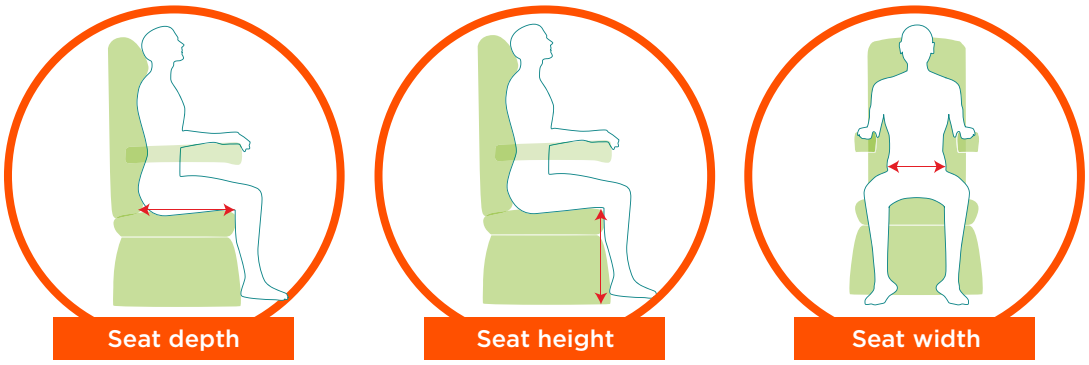


| Chair Size | Part Number | Seat Depth | Seat Height | Seat Width | Weight Capacity |
|---|---|---|---|---|---|
| ☐ SMALL | ☐ CR5404 - Duratek<br>☐ CR5422 - Beige Vinyl<br>☐ CR5443 - Black Vinyl | ☐ 356 mm (14”) #1<br>☐ 406 mm (16”) #2<br>☐ 457 mm (18”)<br>☐ 508 mm (20”) | ☐ 388 mm (15”) #5<br>☐ 406 mm (16”)<br>☐ 432mm (17”) #3<br>☐ 457 mm (18”) #4 | 457 mm (18”) | Seat: 130kg<br>Legrest: 60kg |
| ☐ MEDIUM | ☐ CR5406 - Duratek<br>☐ CR5424 - Beige Vinyl<br>☐ CR5444 - Black Vinyl | ☐ 356 mm (14”) #1<br>☐ 406 mm (16”) #2<br>☐ 457 mm (18”)<br>☐ 508 mm (20”) | ☐ 432mm (17”) #5<br>☐ 457 mm (18”)<br>☐ 482 mm (19”) #3<br>☐ 508 mm (20”) #4 | 508 mm (20”) | Seat: 160kg<br>Legrest: 60kg |
| ☐ LARGE | ☐ CR5408 - Duratek<br>☐ CR5426 - Beige Vinyl<br>☐ CR5445 - Black Vinyl | ☐ 356 mm (14”) #1<br>☐ 406 mm (16”) #2<br>☐ 457 mm (18”)<br>☐ 508 mm (20”) | ☐ 432mm (17”) #5<br>☐ 457 mm (18”)<br>☐ 482 mm (19”) #3<br>☐ 508 mm (20”) #4 | 559 mm (22”) | Seat: 160kg<br>Legrest: 60kg |

| #1 Depth Adjustment ☐ | #2 Depth Adjustment ☐ | #3 Height Adjustment ☐ | #4 Height Adjustment ☐ | #5 Height Reduction ☐ |
|---|---|---|---|---|
| 2 x 2” Kits (100mm) CA2417  | 2” Kit (50mm) CA2416  | 1” Raising Kit (25mm) CA2421  | 2” Raising Kit (50mm) CA2422  | 1” Lowering Kit (18mm) CA2423  |

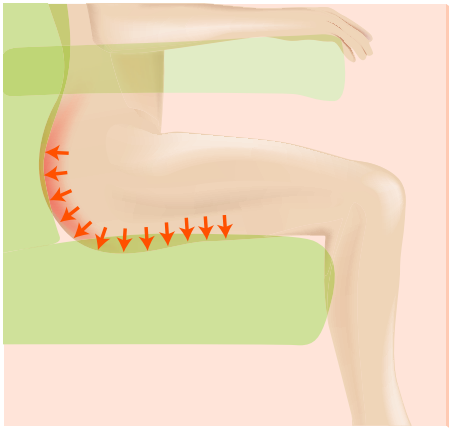


# PRESSURE

**Key thoughts:**

- Selection of pressure cushions is a clinical decision which involves many variables including a risk assessment tool. Professional input from a Tissue Viability Nurse or District Nurse may be required at this stage.
- Vapour permeable fabric as standard on seat base, legrest, backrest and armrests.

## PRESSURE CARE OPTIONS

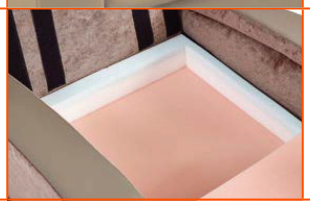
| Description | Long Description | Part Number | | Select |
|---|---|---|---|---|
| Visco (std) | **Prevention:** Heat sensitive memory foam moulds to the user providing maximum support and pressure reduction. Suitable for users who are at risk of pressure damage. 160kg (25 stone) weight capacity. | N/A – standard part |  | ☐ |
| Elata Alternating Air Cushion - E600C | **Management:** powered air cushion for users who have developed pressure damage. Weight capacity: 160 kg. | Fits all sizes E600C |  | ☐ |
| Integrated Self Supply | Allows integration of self-supply cushion of approximate dimensions 17.5" x 17.5" x 3". (Professional input required. See page 11 for instructions). | N/A – standard part |  | ☐ |


## Cushionair:




The Cushionair cushion is a dynamic pressure relieving cushion which is normally recommended for users who have developed pressure damage*.

The principle is that the cells of the cushion inflate and deflate alternately. This provides pressure redistribution via cyclic changes in loading and unloading (http://www.npuap.org/wpcontent/uploads/2012/03/NPUAP_S3I_TD.pdf).

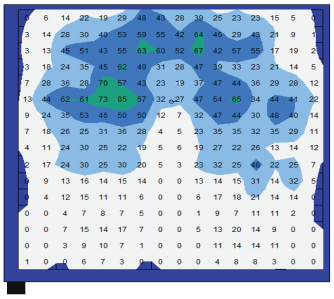


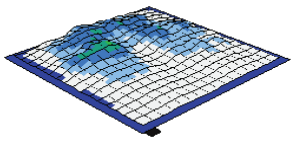

| | |
|---|---|
| Minimum (mmHg) | 0.00 |
| Maximum (mmHg) | 84.57 |
| Average (mmHg) | 21.51 |
| Variance (mmHg²) | 347.67 |
| Standard deviation (mmHg) | 18.65 |
| Coefficient of variation (%) | 86.69 |
| Horizontal center (cm) | 21.73 |
| Vertical center (cm) | 28.02 |
| Sensing area (cm²) | 1866.24 |
| Regional distribution (%) | 100.00 |

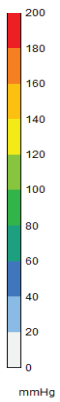


*This is based around anecdotal evidence and Accora's experience; we do not yet have clinical studies to evidence this although these are in progress.

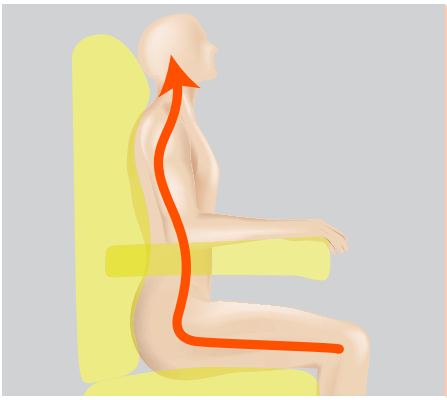

# POSTURE

**Key thoughts:**

- Tilt-in-space can significantly improve midline sitting so increased lateral support may not be necessary if the user tends to slump mildly to one side.
- The Multi-adjustable Pillow Backrest is more adjustable than the Lateral Support Backrest and Adjustable Lateral Support System because filling can be removed to accommodate different postural requirements.

## POSTURE MANAGEMENT OPTIONS

| Description | Long Description | Part Number | | Select |
|---|---|---|---|---|
| Multi-adjustable Pillow Backrest (std) | **Accommodation:** accommodates posture that has already developed i.e. kyphosis; our most popular backrest. | N/A Standard Part |  | ☐ |
| Lateral Supports | **Accommodation:** used under the Multi-adjustable Pillow Backrest cushions to provide moderatepostural support **SOLD IN A PAIR** | CA2413 |  | QTY 2 ☐ |
| Cocoon Backrest | **Support:** offering additional comfort and support, adjusted by changing the position or amount of stuffing in any of its 9 cells. The "reduced shear backrest" mechanism on the Comfort chair complements this backrest by providing extra trunk support, whilst retaining comfort and pelvic position. | Fits all sizes **Beige** - CA2713 |  | ☐ |
| | | **Black** - CA2711 |  | ☐ |
| Lateral Support Backrest | **Encouragement:** encourages a good midline sitting position for users whose posture is flexible. | **Beige Small** - CA2406<br>**Beige Medium** - CA2407<br>**Beige Large** - CA2408 |  | S ☐<br>M ☐<br>L ☐ |
| | | **Small** - CA2540<br>**Medium** - CA2541<br>**Large** - CA2542 |  | S ☐<br>M ☐<br>L ☐ |
| Adjustable Lateral Support Backrest | **Management/correction:** contoured foam backrest with adjustable Velcro Lateral Supports. Provides more aggressive lateral support and mild positioning/ correction. Supplied with breathable, pressure reducing fabric. | **Small** - CA2410<br>**Medium** - CA2411<br>**Large** - CA2412 |  | S ☐<br>M ☐<br>L ☐ |

continued overpage




| | | | | |
|---|---|---|---|---|
| Profiled Headrest | Suitable for users with reduced lateral head control. Fits over the shoulders. | Fits all sizes<br>**Beige** - CA2419 |  | ☐ |
| | | **Black** - CA2530 |  | ☐ |
| 2-Button Handset | Suitable for Cognitive Impaired person e.g. dimentia (Chair comes Standard w/ 4-button Handset). | Fits all Chair Sizes<br>SC2316 |  | ☐ |
| Legrest Channel | Provides side support on legrest to reduce abduction | **Small** - CA2520<br>**Medium** - CA2521<br>**Large** - CA2522 |  | S ☐<br>M ☐<br>L ☐ |

| Description | Long description | Part Number | | Select |
|---|---|---|---|---|
| Right-hand Side Configura Handset Holder | A great addition for users with dementia, or similar conditions, who may have difficulties locating the handset. Left or right-hand side option. Ideal for users with limited movement on one side of their body from conditions like brain injury or stroke. Users can easily reach the handset from any position, including recline and tilt-in-space, while maintaining optimal posture. | CA2556 |  | ☐ |
| Left-hand Side Configura Handset Holder | | CA2554 |  | ☐ |

## TRANSFER OPTIONS

| Description | Long description | Part Number | | Select |
|---|---|---|---|---|
| Lefthand side Dropdown Armrest | **Transfer**<br>Enables user to do a slideboard transfer in and out of the chair with ease. Can be bolted on at any time if the users needs change. | **Small**<br>**Duratek** - CA2484<br>**Beige Vinyl** - CA2492<br>**Black Vinyl** - CA2534<br>**Medium / Large**<br>**Duratek** - CA2488<br>**Beige Vinyl** - CA2496<br>**Black Vinyl**- CA2535 |  | S ☐<br>M ☐<br>L ☐ |
| Righthand side Dropdown Armrest | **Transfer**<br>Enables user to do a slideboard transfer in and out of the chair with ease. Can be bolted on at any time if the users needs change. | **Small**<br>**Duratek** - CA2486<br>**Beige Vinyl** - CA2494<br>**Black Vinyl** - CA2536<br>**Medium / Large**<br>**Duratek** - CA2490<br>**Beige Vinyl** - CA2498<br>**Black Viny** - CA2537 |  | S ☐<br>M ☐<br>L ☐ |

# EXTRAS
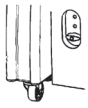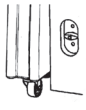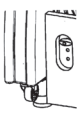
| Description | Long description | Part Number | | Select |
|---|---|---|---|---|
| Armrest overhang on outside (std) | Provides a good grip to assist the user rolling from side to side for easier fitting of slings. | Standard Part |  | ☐ |
| Armrest overhang on inside | Reduces width between armrests by approximately 38 mm (1.5”). Provides extra arm and elbow support if the user is ‘pear shaped’ or very thin. | Standard Part |  | ☐ |
| Right hand control (std) | N/A | Standard Part |  | ☐ |
| Left hand control | N/A | Standard Part |  | ☐ |
| Carer control | Fitted to backrest out of user reach | Standard Part |  | ☐ |
| Rise (std) | 23° of tilt. The most frequently used rise setting. This provides a gentle, progressive movement. | Standard Part |  | ☐ |
| Medium rise | 14° of tilt. There is not a specific user group for this setting – it is normally used if ‘Rise’ or ‘Flat Rise’ is inappropriate. | Standard Part |  | ☐ |
| Flat rise | 2° of tilt and 254 mm (10”) of rise, used for users with Muscular Dystrophy who need to lock knees before transferring. (Not available on Small chair) | Standard Part |  | ☐ |
| Legrest Setting 1. | Standard Position | Standard Part |  | ☐ |
| Legrest Setting 2. | Standard Position + 25 mm (1”) | Standard Part |  | ☐ |
| Legrest Setting 3. | Standard Position + 50 mm (2”) | Standard Part |  | ☐ |
| Legrest Setting 4. | Standard Position -25mm (1”) | Standard Part (rotate the footboard 180 degrees) |  | ☐ |

# TECHNICAL SPECIFICATION

- To provide comfortable seating
- To raise a user from a sitting to a standing position
- To lower a user from a standing to a sitting position
- To tilt in space a user

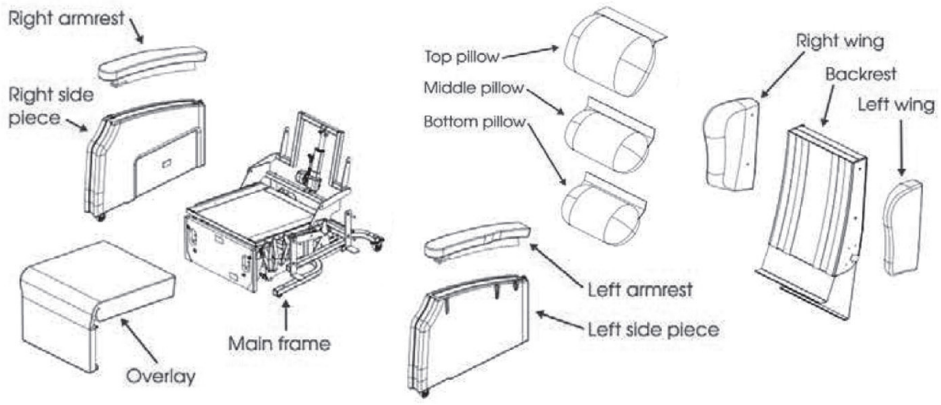



Figure 1. Spare parts of Configura® Comfort 1 Chair

## OVERALL DIMENSIONS

| | SMALL | MEDIUM | LARGE |
|---|---|---|---|
| WIDTH | 711 mm (28”) | 762 mm (30”) | 812 mm (32”) |
| LENGTH - UPRIGHT POSITION | 940 mm (37”) | 940 mm (37”) | 940 mm (37”) |
| LENGTH - RECLINE POSITION | 1702 mm (67”) | 1702 mm (67”) | 1702 mm (67”) |
| HEIGHT | 1092 mm (43”) | 1143 mm (45”) | 1143 mm (45”) |
| MASS | 65 kg (143 lbs) | 68 kg (150 lbs) | 71 kg (156 lbs) |

## SPECIFIC DIMENSIONS

| | SMALL | MEDIUM | LARGE |
|---|---|---|---|
| SEAT HEIGHT | 406 mm (16”) | 457 mm (18”) | 457 mm (18”) |
| SEAT HEIGHT WITH (HAKS) | 457 mm (18”) | 508 mm (20”) | 508 mm (20”) |
| SEAT WIDTH | 457 mm (18”) | 508 mm (20”) | 558 mm (22”) |
| SEAT DEPTH POSITION | 457 mm (18”) | 508 mm (20”) | 508 mm (20”) |
| ARMREST HEIGHT | 180 mm (7”) | 180 mm (7”) | 180 mm (7”) |
| RISE HEIGHT TO FRONT OF CHAIR | | | |
| STANDARD | 546 mm (21.5”) | 508 mm (20”) | 508 mm (20”) |
| MEDIUM | 622 mm (24.5”) | 572 mm (22.5”) | 572 mm (22.5”) |
| FLAT | NA | 711 mm (28”) | 711 mm (28”) |
| CHAIR TO WALL / FURNITURE | 457 mm (18”) | 483 mm (19”) | 483 mm (19”) |
| SAFE WORKING LOAD | 130 kg | 160 kg | 160kg |
| BACK HEIGHT | 700 mm | 700 mm | 700 mm |

# CLEANING INSTRUCTIONS FOR THE FABRIC

Clean and care for the fabric on a regular basis and it will maintain its condition and appearance.

The new generation of Microfibre cleaning cloths without soap are extremely effective for general care and maintenance. You can also wipe clean regularly with a damp cloth to remove dust particles and light soiling. It is recommended to vacuum the upholstery once a week.

Most stains can be removed by wiping clean. Remove spillages promptly (within 15 minutes) with cold water. For extreme soiling control use a 5% Sodium Hypochlorite (Bleach) solution.

Do not machine wash or dry clean.

Do not apply other chemicals such as polish, solvents or stain remover.

# INTEGRATED SELF-SUPPLY OVERLAY

The integrated Self Supply Overlay comes as standard on all chairs.
This allows for easy integration of existing cushions.
The size of the aperture is 445 mm x 445 mm x 76 mm (17.5” x 17.5” x 3”).

This means the height of the seat and armrests remain unchanged when using another cushion.



## STEP 1

Remove the backrest cushions.



## STEP 2

Unzip the cushion cover to expose the foam.



## STEP 3

Pull the foam out carefully, starting at the rear and cut the perforations at the front.



## STEP 4

This leaves the cut out for the existing cushions to be inserted. Insert the cushion and re zip the cover.

At Enable we source quality, affordable and innovative products globally and have them available for fast delivery. But we do more than just provide products.

We give the knowledge and support to help our customers grow their business and deliver better life care.

enable lifecare
Fast intelligent care

ENABLE LIFECARE PTY LTD
T 1300 370 370
enablelifecare.com.au